

ADVANCE Bath Chair

LECKEY



Provides the user with the level of postural support necessary to ensure their safety and protection during bathing and outdoor leisure. It is configurable in several positions depending on the postural needs of the user and adjustable in height to adapt to the needs of the family member or caregiver.

The height adjustment of the legs allows you to configure the hammock at the appropriate height to work in a comfortable and ergonomic position, reducing the risk of strain or back injuries. **Image 1.**

It allows tilting, and the lower extremity support can be adjusted in 4 positions. **Image 2.**

Optional wheel base. **Image 3.**

Light and robust chassis, made of stainless aluminum with steel and plastic elements. **Image 4.**



Image 1



Image 2



Image 3

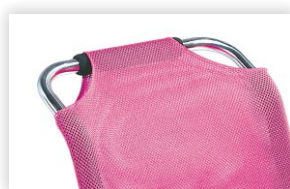


Image 4

Keeping you in mind. More options.

Equipment manufactured and adapted for children and adults with disabilities.

SPECIFICATIONS

Sizes	4
Aproximated age	1 to 18 years
Max user weight	72 kg / 158 lbs
User height	75 to 165 cm / 29.5" to 65"
Backrest recline	0° - 90°
Accessories	Optional
Base width	60cm/23.6"
Total base length	90 to 115 cm / 35.4" to 45.3"
Fabric color	Blue and pink

Required user
measurement
(Total height)



BENEFITS

Suitable for patients with neuronal and psychomotor damage or cognitive psychic disability.

Discover our products



DESIGNED FOR OUR TERRITORIES

Our products are designed to meet the needs of each of our territories, offering quality and affordable prices.



FUNCTIONAL VERSALITY

We have a variety of products ranging from a wide range of manual and power wheelchairs and mobility scooter.

**For more information on obtaining this product
Contact your local Loh Medical representative**

Find your local representative at www.lohmedical.com   

ver_2025.11.11E

