

Stander Prone & Supine



Prone

Supine



Image 1



Image 2



Image 3

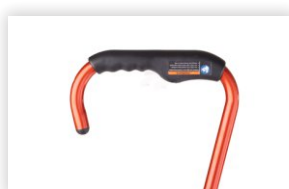


Image 4

Multi-position stabilizer that can easily switch between face down and face up, giving you the most flexibility and the most options for the user. Provides support and stability while encouraging weight bearing and use of hands and arms.

It can be tilted 90° from horizontal to completely vertical. Adjust the ideal angle for each user. **Image 1.**

Body support with large, padded area, which supports the user's pelvis and trunk with optional Split Knee Support. **Image 2.**


Independent hip abduction with 0-30° of independent abduction for each lower extremity. **Image 3.**





The caregiver handle allows for easy maneuvering and can be quickly adjusted up or down. **Image 4.**



Keeping you in mind. **More options.**

Equipment manufactured and adapted for children and adults with disabilities.

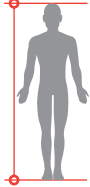SPECIFICATIONS	
ITEM	S1 S410 S2 S420 S430
Size	S1, S2 & S3
Frame	Aeronautical aluminum
User measurement	28" a 65" / 71 a 165 cm
Max user weight	185 lbs / 84 kg (S3)
Product weight	Prone 47, 72 & 104 lbs / 22, 33 & 47 kg Supine 52, 78 & 115 lbs / 24, 35 & 52 kg
Accessories	Optional
Settings	No tools required
Frame colors	

Required user measurement (Total height)	
	Option for growth
Adjustable sandals Rotate 180°	
	Foot brakes on both side wheels

BENEFITS

Allows upright position for adults and children with physical disabilities.

Discover our products



DESIGNED FOR OUR TERRITORIES



Our products are designed to meet the needs of each of our territories, offering quality and affordable prices.



FUNCTIONAL VERSALITY

We have a variety of products ranging from a wide range of manual and power wheelchairs and mobility scooter.

**For more information on obtaining this product
Contact your local Loh Medical representative**

Find your local representative at www.lohmedical.com   

ver_2025.11.11E

